


STEAM INTEGRATOR X2

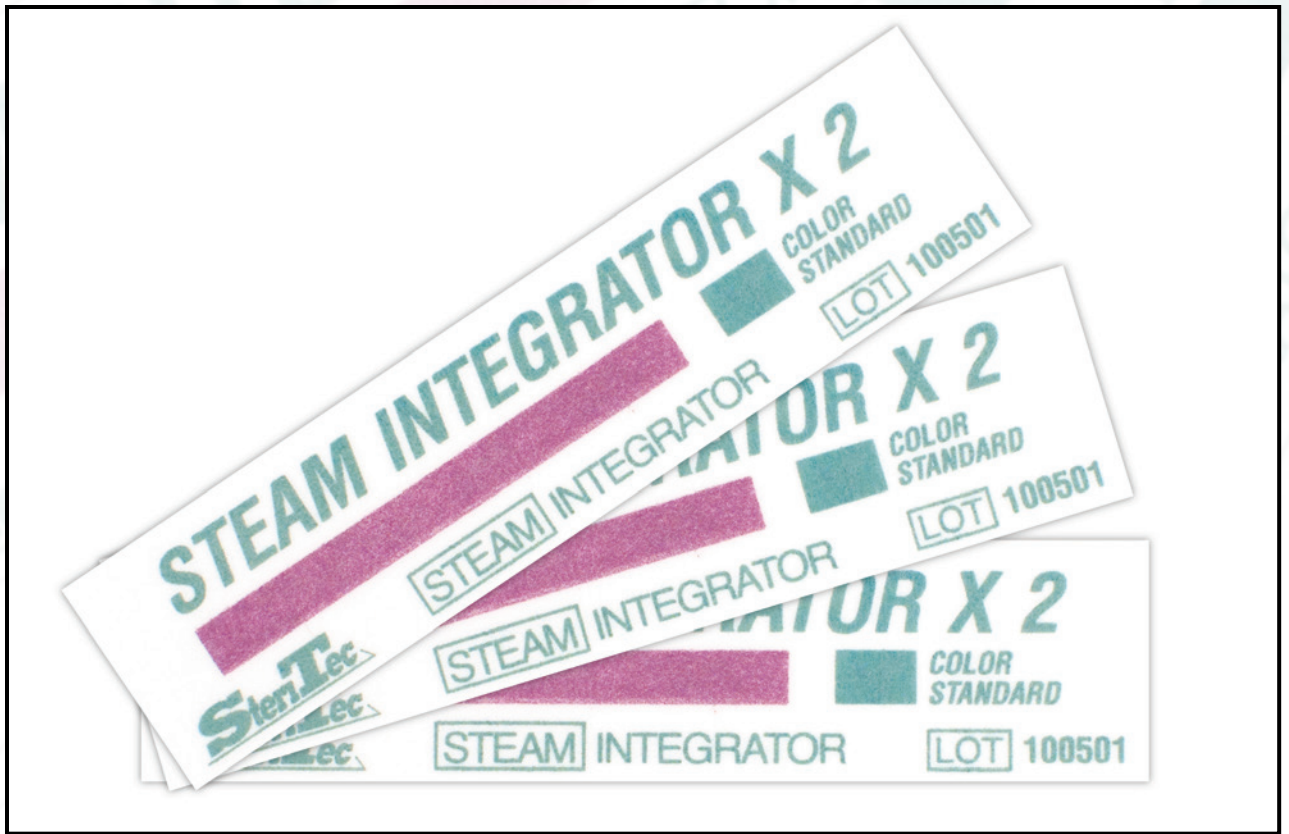

Steam Integrator X2

Chemical Integrator for Steam Sterilization

Catalog No. CI 143

- Vivid purple-to-green color change
 - Full cycle monitor that is sterilant specific meaning it must have steam to change
 - Will not complete a color change if excess air is present
 - Turns dark brown under superheated steam conditions
 - Designed for steam sterilization at 132°C (270°F)
 - Color change is designed to occur after twice biological kill has been met
- Released for market by the FDA as an Integrator
- Exceeds performance requirements of integrators described by both ANSI/AAMI and EN ISO11140-1
- Lead free

ORDERING INFORMATION

PART #	PACKAGED
CI 143	1 CASE = 20 BAGS OF 250

SteriTec
www.steritecproducts.com